


# Wallace® Amniocentesis Needles



Spiral engraving for enhanced visualization

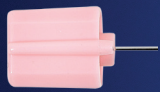

Stylet to prevent blockage



Amniocentesis needle 9022 (90 mm, 22g)



Amniocentesis needle 12021 (120 mm, 21g)



Amniocentesis needle 15021 (150 mm, 21g)

This is a high quality needle specifically designed for optimum handling and confidence of tip placement during amniotic fluid sampling.

The single-use amniocentesis needles are provided in individual sterile packs for the application of amniotic fluid sampling. The needle tips have spiral engraving to provide clear visibility under ultrasound guidance. Each needle is provided with a self locating stylet to ensure the lumen is free of blockages once the tip is in place.

## Exceptional Performance, Safety and Control:

- 1cm of spiral engraving at needle tip provides a reflective surface for enhanced visualization under ultrasound guidance
- Sharp, non-coring needle tip to minimize patient trauma
- Self-locating stylets to prevent blockage
- Color coded hubs for easy identification of needle length and gauge
- Range of needle lengths: 90mm, 120mm, 150mm

Product Code	Description
9022	Amniocentesis needle 90mm, 22g – 10/Box
12021	Amniocentesis needle 120mm, 21g – 10/Box
15021	Amniocentesis needle 150mm, 21g – 10/Box



To place an order, or to learn more, contact your CooperSurgical representative, visit [www.CooperSurgical.com](http://www.CooperSurgical.com), or call 800.243.2974 or 203.601.5200.